


Cleaning the Sample Compartment

Waters
THE SCIENCE OF WHAT'S POSSIBLE.™

- This section covers:
 - Removing the front panel
 - Removing the carrier
 - Replacing the carrier and front panel


Waters
THE SCIENCE OF WHAT'S POSSIBLE.™

Remove the Carousels

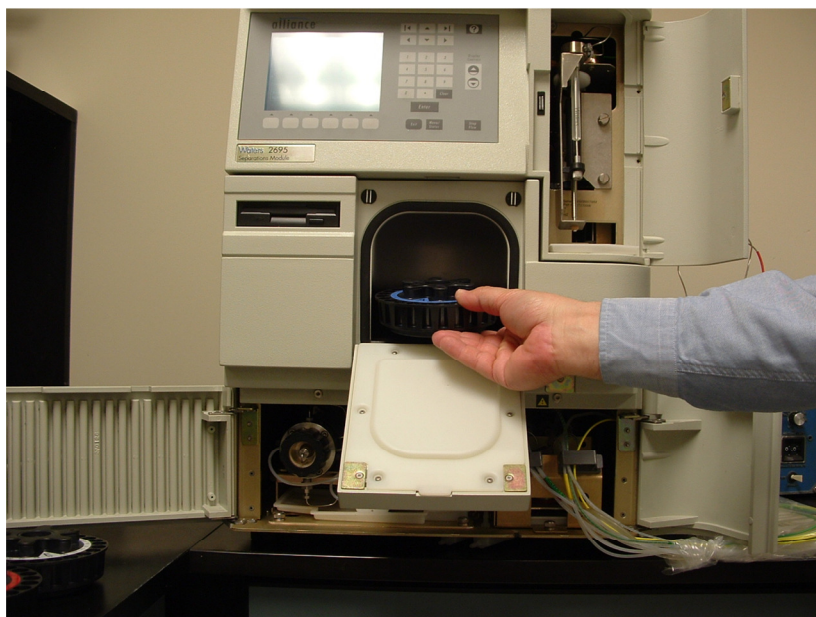


- Open the solvent conditioning, solvent delivery, and syringe compartment doors.
- Open the sample compartment door and remove the carousels from the compartment using the carousel advance function on the display window.
- After removing all the carousels, bring the carrier to the "A" position

Removing the Front Panel

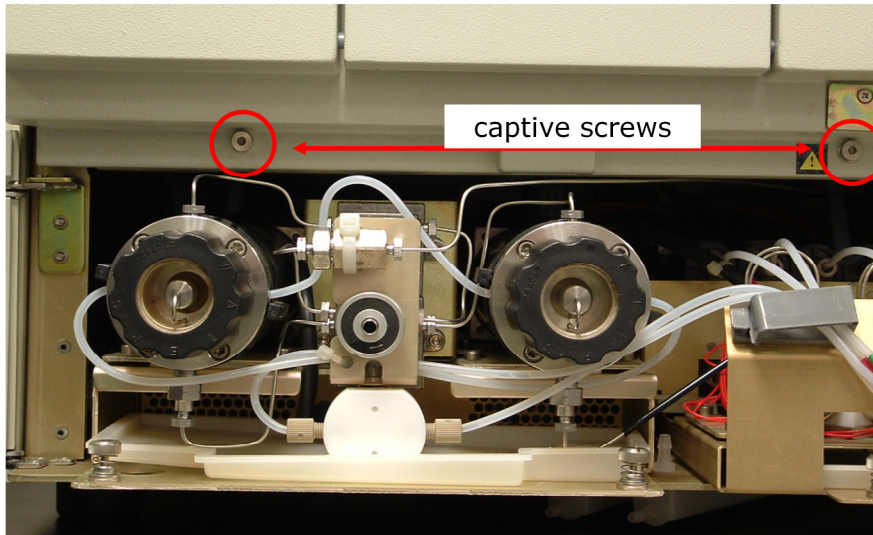

- There are two captive Allen screws (see above) that hold the front panel in place.

Loosen the Captive Screws

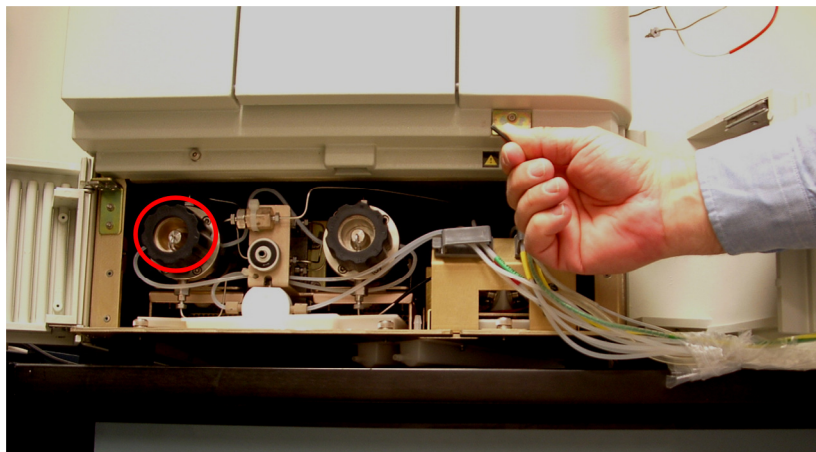


- Loosen the captive screws using a 3 mm Allen wrench.

Remove the Front Panel




- Pull the bottom of the front panel forward and lower it away from the instrument

Remove Carrier

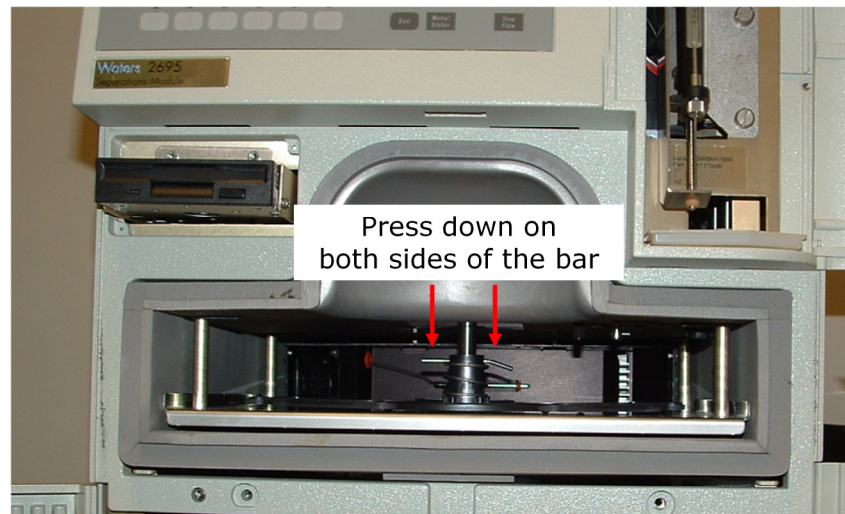


- Press down on both sides of the metal bar as indicated above.
- Slide the carrier towards you.

Clean Sample Compartment

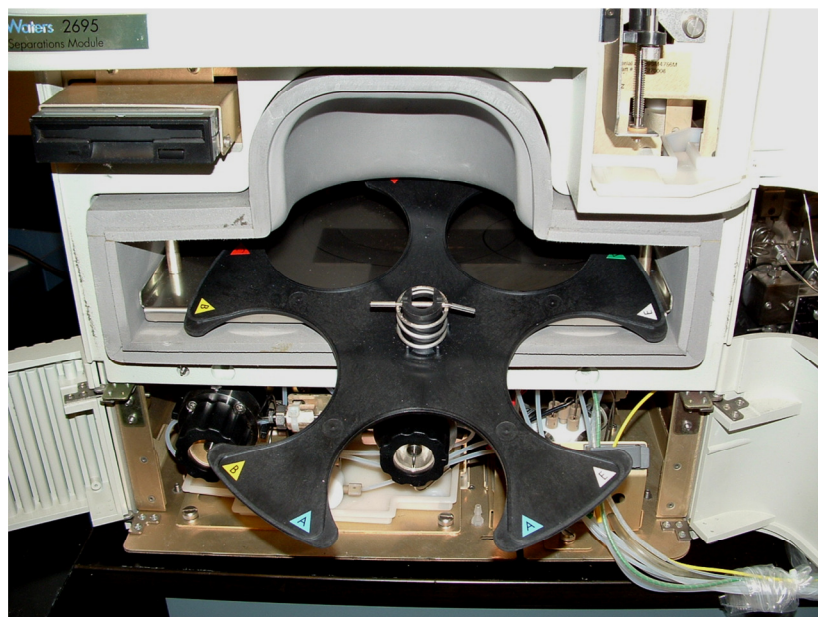

- Once the carrier is removed, the carrier and compartment can be cleaned using an appropriate cleaning solvent.
- After cleaning the sample compartment and carrier, replace the carrier.
- The "A" carousel position must be facing the front of the compartment.

Replace the Carrier

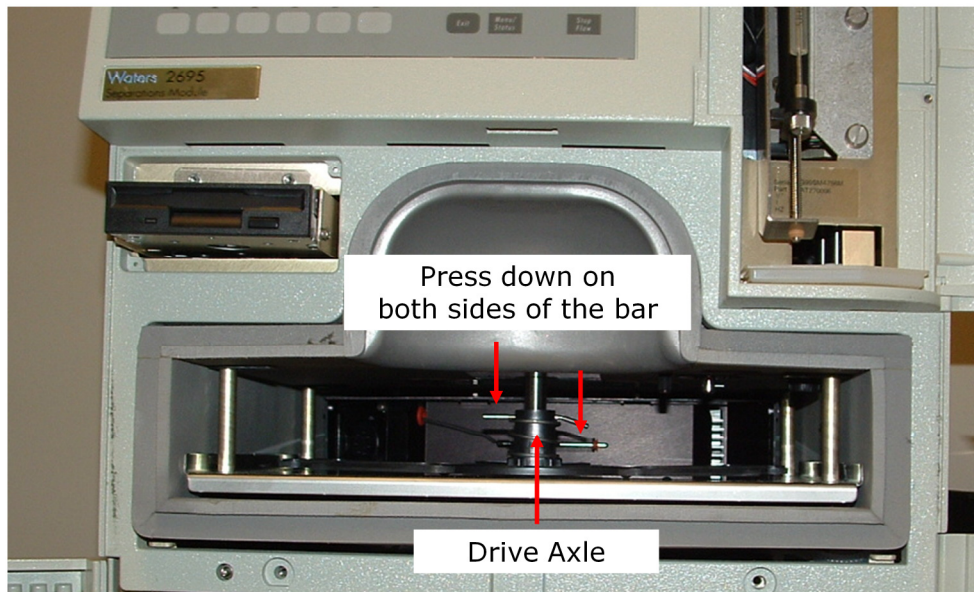

- Press down on both sides of the metal bar indicated above and slide the carrier over the carrier drive axel.


Reinstall the Front Panel



- Place the top of the front panel in place and swing the bottom in towards the instrument.

Note: It may require some wiggling to get the front panel over the diskette drive.

Tighten the Captive Screws

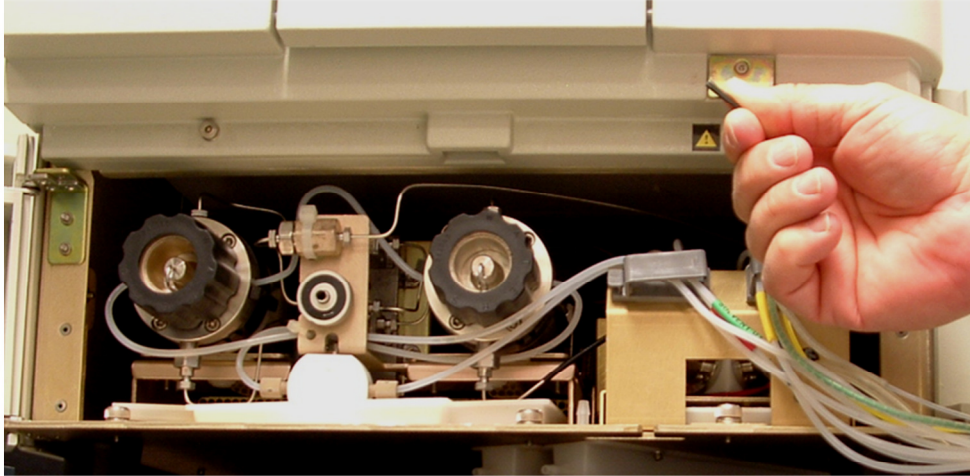

- Tighten both captive 3 mm Allen screws to secure the front panel to the instrument.
- Replace the carousels.